GB

INSTRUCTION MANUAL

Bathroom Scale

Model: PW-304

Thank you for choosing Vox Electronics bathroom scale. To ensure continued operation and accuracy, please read these operating instructions carefully before use. By passing certain procedures may result in incorrect readings.

IMPORTANT SAFETY PRECAUTIONS


1. Never place the scales on carpet or on uneven floor.
2. The scales are not waterproof; never immerse the scales in water. To clean the surface, use a damp cloth or a glass cleaner. Do not use soap or other chemicals. Please keep the scales away from water, heat and extreme cold.
3. The scales are a high precision measuring device. Never jump or stamp on the scales or disassemble. Please handle the scales carefully to avoid breakage.
4. Remove the batteries if the scales are not being used for a long period of time.
5. Please check the batteries if the scale malfunction. Change to new batteries if required.
6. Please remove insulating strip on the batteries before use.
7. If batteries leak, change them immediately they may damage the scales.
8. To ensure weighing accuracy, when the scales are first used please place them on an even surface and wait for 10 seconds before stepping on the scales. This also applies if the scales are moved.
9. To ensure accuracy, we suggest you take your measurements at the same time every day.

SPECIFICATION

1. Auto ON
2. Auto OFF
3. High precision strain gauge sensors system
4. Weight Units: kg, lb, st: lb
5. Maximum Capacity: 150kg/330lb/24st
6. Tempered glass platform: 312x312x6mm
7. LED size: 86x41mm
8. Graduation: 0.1kg/0.2lb/0.2lb(1lb \geq 20st)
9. Overload indication: ERR
10. Low batteries indication: LO
11. Power supply: 4 \times AAA batteries
12. Batteries included
13. 1 year manufacturer's warranty

OPERATING INSTRUCTIONS

1. Remove all packaging before use.
2. Locate the 4 \times AAA batteries and the battery holder at the back of the scales, open the cover, install the batteries into the holder and close the cover. Please ensure that the batteries are placed the correct way as indicated in the holder.

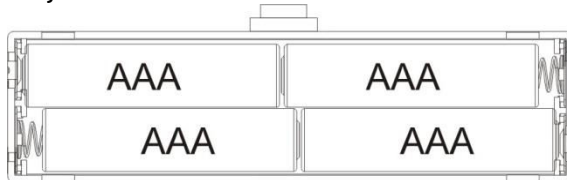3. Press UNIT key to switch on and select the unit measurement.
4. Place the scales on solid, flat ground, let the scales turn off automatically.
5. Step on the scales gently, position your feet evenly near the center and keep your body still and erect until the scales stabilize your weight. The scales will flick and lock on the LED screen after your weight is steady, then automatically turns off after 10 seconds.
6. The scales will automatically turn off if no further use required.

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