80mm

80mm

80mm

125mm

GBR

ITA

DEU

RUS

SRB

ESP

PRT

GRC

SVN

HRV

BIH

MNE

MKD


ALB

**PW 304**


OPERATING INSTRUCTIONS BATHROOM SCALE  
MANUALE DI ISTRUZIONI BILANCIA PESAPERSONE  
BEDIENUNGSANLEITUNG PERSONENWAAGE  
ИНСТРУКЦИЯ ПО ЭКСПЛУАТАЦИИ **ВЫГОВЫЕ  
НАПОЛЬНЫЕ ВЕСЫ**  
UPUTSTVO ZA UPOTREBU TELESNA VAGA  
MANUAL DE USUARIO PESO CORPORAL  
MANUAL DO USUÁRIO BALANÇA DIGITAL  
ΟΔΗΓΙΕΣ ΧΡΗΣΗ ΖΥΓΑΡΙΑ  
NAVODILA ZA UPORABU OSEBNA TEHTNICA  
URUTE ZA UPORABU OSOBNA VAGA  
UPUTSTVO ZA UPOTREBU TJELESNA VAGA  
УПАТСТВА ЗА РАКУВАЊЕ **ТЕЛЕСНА ВАГА**  
MANUALI I PERDORUESIT PESHORE TRUPORE

80mm                      80mm                      80mm

**Function keys**

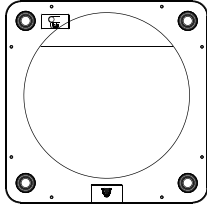1. Remove all packaging before use.
2. Locate the AAA battery and the battery holder at the back of scale, open the cover, install the batteries into the holder and close the cover.
3. Please ensure that batteries are placed in the correct way as indicated in the holder

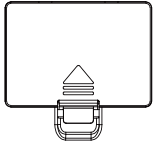4. Press UNIT key to switch on and select unit measurement.
5. Place the scales on solid, flat ground, let the scale turn off automatically.
6. Step on the scale gently, position your feet evenly near the center and keep your body still and erect until the scale stabilize your weight. The scale will flick and lock on the LCD

screen after your weight is steady, then automatically turn off after 10 seconds.  
7. The scale will automatically turn off if no further use required.

**Error indicators**

**Overload:** 


**Low battery supply:** 

**Working conditions:**

1. Operating temperature range: 0 ~ 40 °C
2. Operating humidity range: ≤ 90% RH
3. Storage temperature: -20 °C ~ 60 °C

**Guarantee**

These scales have been checked for accuracy prior to leaving the factory and have a guarantee against any defect in either workmanship or materials for a period of one year.

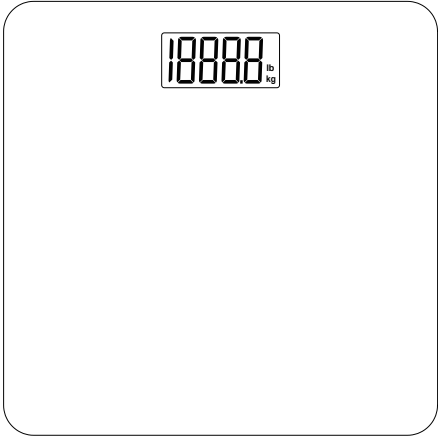**Bathroom Scale  
User Manual**

PW 304


125mm

To ensure continued operation and accuracy, please read these operating instructions carefully before use. By passing uncertain procedures may result in incorrect readings. We hope that you enjoy using your new scales.


**Important safety precautions**

When using electrical appliances, certain basic safety precautions must be followed, including the following:

1. This product must not be used by children or by other people without assistant or supervision if their physical, sensory or mental condition limits safe use of the appliance.
2. In order to ensure reliable readings and guarantee connectivity, please stand on the scales with bare feet. Please make sure that both the scales and your feet are completely dry before use.
3. Never place the scales on carpet or on uneven floor.
4. The scales are not waterproof, never immerse the scales in water. To clean the surface, use a damp cloth or a glass cleaner. Do not use soap or other chemical. Please keep the scales away from water, heat and extreme cold.
5. The scales are a high precision measuring device. Never jump or stamp on the scales or disassemble. Please handle the scales carefully to avoid breakage.
6. The scales are for domestic use only. You should consult a doctor when undertaking any diet or exercise program.
7. Remove the batteries if the scales are not being used for a long period of time.
8. Please check the batteries if the scale malfunctions. Change to new batteries if required.
9. Please remove insulating strip on the batteries before use.
10. If batteries leak, change them immediately as they may damage the scales.

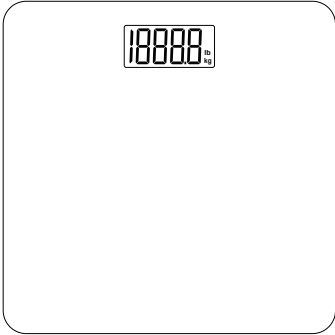**Full display**

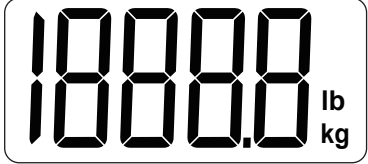
**Features & Specification**

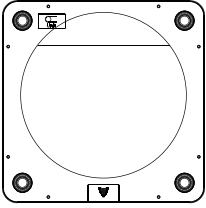1. Auto on
2. Auto off
3. High precision strain gauge sensors system
4. Weight Units: kg, lb
5. Maximum Capacity: 180kg / 396lb
6. Material: Tempered glass platform
7. Product size : 300 x 300 x 23.5mm
8. LCD / LED display
9. Graduation: 0.05kg / 0.1lb
10. With room temperature
11. Overload indication: "Err"
12. Low batteries indication: "Lo"
13. Power supply: 2 or 3 x AAA batteries ( optional USB charge battery )
14. Batteries included
15. 1 year warranty

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